

Oro-nasal masks for hospital use



The solution
to benefit from a
complete expertise
in non-invasive
ventilation

The Respireo™ Hospital F oro-nasal masks represent a range of **disposable** and **reusable** masks, **vented** and **non vented**, specifically dedicated to **hospital use** and entirely designed for **non-invasive ventilation**.

Available in **three different sizes**, these masks are designed to provide a comfortable fit for a wide range of patients and thus maximizing treatment tolerance.

Respireo Hospital F

A simple, thoughtful and safe design

Adjustable headgear

Designed to be easy to set-up and to fit all head morphologies

Silicone double-wall cushion

Provides a comfortable and optimal seal by adapting to the patient's face

Quick-release clips

Designed for easy headgear attachment and quick release in case of emergency

Silicone forehead pad

Designed to ensure mask stability and patient comfort

Adjustment knob

Designed to adjust the pad on the patient's forehead (16 possible positions)

Transparent frame


Designed to provide an optimal line of sight for the clinicians

Pressure port

Single port for measuring pressure

360° rotating elbow

Facilitates the movement and the position of the tube

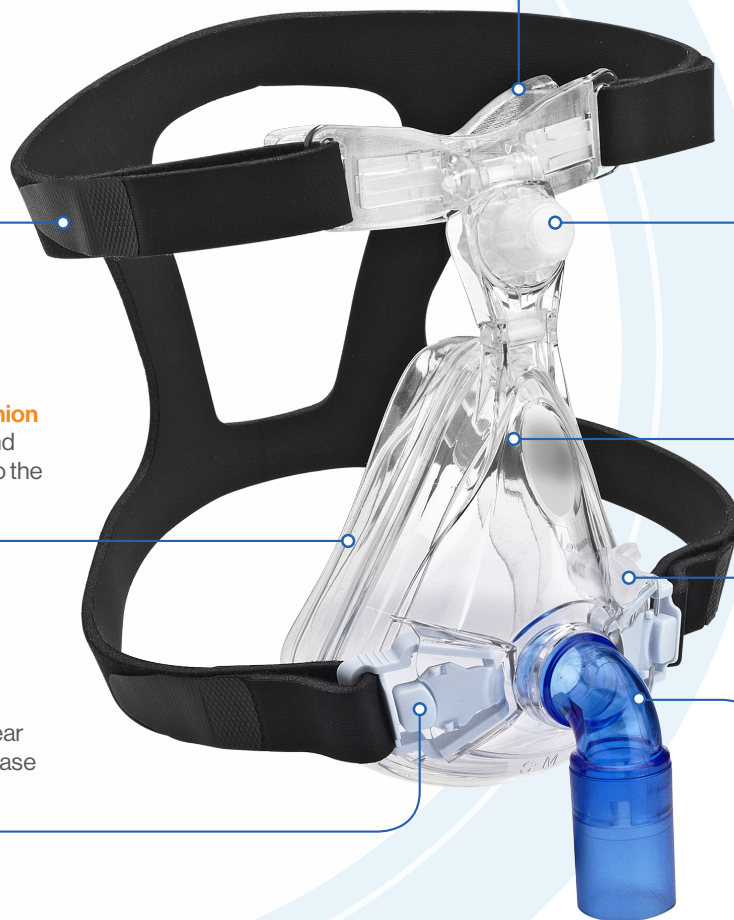

A wide range of masks for a variety of patients

Disposable
single use for up to 7 days



Non vented with female connector



Non vented with anti-asphyxia valve and female connector



Vented with male connector



- 3 cushion sizes available (S, M, L)
- Reclining forehead support for an optimal stability
- Package including mask, cushion, headgear and sizing gauge

Reusable
multi-patient usedays




Non vented with female connector


Fitting instructions for Hospital F masks

First, use the sizing-gauge on the package to determine the correct mask size for the patient


1 Turn the knob to the maximum setting

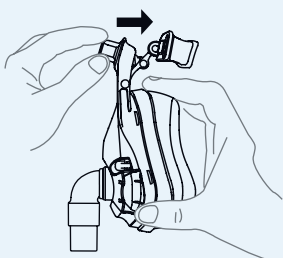


2 Place the mask on the patient's face

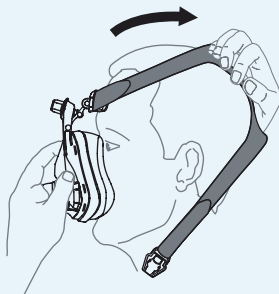

3 Smoothly adjust the headgear

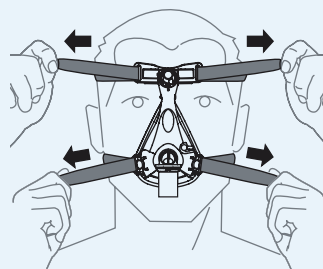

4 Adjust the knob to obtain a better seal

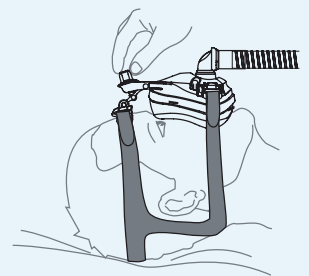

Product references

Respireo™ Hospital F Disposable

Product Reference	Description	Size	Distributor code
472539	Respireo™ Hospital F Vented Disposable	S	KM226400
472540	Respireo™ Hospital F Vented Disposable	M	KM226500
472541	Respireo™ Hospital F Vented Disposable	L	KM226600
472542	Respireo™ Hospital F Non Vented Disposable	S	KM222100
472543	Respireo™ Hospital F Non Vented Disposable	M	KM222200
472544	Respireo™ Hospital F Non Vented Disposable	L	KM222300
472526	Respireo™ Hospital F Non Vented Disposable with anti-asphyxia valve	S	KM229132
472527	Respireo™ Hospital F Non Vented Disposable with anti-asphyxia valve	M	KM229133
472528	Respireo™ Hospital F Non Vented Disposable with anti-asphyxia valve	L	KM229134

Respireo™ Hospital F Reusable

Product Reference	Description	Size	Distributor code
472545	Respireo™ Hospital F Non Vented Reusable	S	KM218800
472546	Respireo™ Hospital F Non Vented Reusable	M	KM218900
472547	Respireo™ Hospital F Non Vented Reusable	L	KM219000

Spare parts:

472453	Hospital headgear x10	One size	KM222400
472615	Cushion Respireo™ Primo F/Hospital F x10	S	KM218500
472616	Cushion Respireo™ Primo F/Hospital F x10	M	KM218600
472617	Cushion Respireo™ Primo F/Hospital F x10	L	KM218700


Manufactured by

AIR LIQUIDE MEDICAL SYSTEMS S.r.l
Via Bisceglie, 66 - 20152 Milano (MI) - ITALY
PLANT: Via dei Prati, 62 - 25073 Bovezzo (BS) - ITALY
Tel. +39 030.201.59.11
Fax +39 030.209.83.29

Distributor Contact

AIR LIQUIDE MEDICAL SYSTEMS S.A.
Parc de Haute Technologie
6, rue Georges Besse
92182 Antony Cedex
France

medicaldevice.airliquide.com

 Hotline ALMS: +33 179 51 70 01



Respireo™ Hospital F Vented disposable is an oronasal mask with calibrated exhalation orifices intended to be used for short time periods (less than 7 days) in a hospital setting by a single adult patient (weighing over 30 kg), who has been prescribed for a non-invasive positive pressure ventilation therapy (NIPPV), e.g. continuous positive airway pressure therapy (CPAP) or Bilevel positive airway pressure therapy (BiPAP).

Respireo™ Hospital F Non Vented disposable is an oronasal mask without calibrated exhalation orifices intended to be used, for short time periods (less than 7 days) in a hospital setting, by a single adult patient (weighing over 30 kg) in controlled non-invasive positive pressure ventilation (NIPPV).

Respireo™ Hospital F Non Vented disposable with anti-asphyxia valve is an oronasal mask without calibrated exhalation orifices and with an anti-suffocation valve, intended to be used for short time periods (less than 7 days) in a hospital setting by a single adult patient (weighing over 30 kg) in controlled non-invasive positive pressure ventilation (NIPPV).

Respireo™ Hospital F Non Vented reusable is a reusable oronasal mask without calibrated exhalation holes intended for use in hospital settings by multiple adult patients (weighing more than 30 kg) under Non-invasive Positive-Pressure Ventilation (NPPV).

Medical device class IIa CE0051. Please read the user manual.